

# Made to measure

## Ordering sheet

for made to measure products (round knitted)



memory<sup>®</sup>

**Elegant**  
medical compression

lastofa<sup>®</sup>

**Universal**  
medical compression

ofa bamberg

we care for you

# The correct method of taking measurements prior to fitting compression stockings

**In order to achieve a good fit and the medically correct amount of pressure on the leg, it is essential that measurements be taken duly and accurately. Stockings that fit well prevent the development of congestion and improve the acceptance of a compression therapy among patients.**

Measurements must be taken on an unclothed leg that is free of oedemas (finger pressure test in the ankle region: an indent will remain in the event of an oedema). The ideal time to take measurements is thus in the morning or subsequent to the removal of the compression bandage.

The oedema-free leg should be measured whilst the patient is in a standing position. The leg that is measured must be the standing leg. If the patient is to be fitted with stockings for both legs, each leg must be measured individually.

The circumference should be measured first, whereby the various measuring points of the circumference measurements should be marked using a pen suitable for use on the skin. The length can then be measured easily using these marked points.

**Measurements should be taken in an upward direction at exactly the leg measuring points stipulated in the table.**

Measuring points	Explanation
$\ell$	Length, measured from the heel to the appropriate measuring point
$c$	Circumference
T	Waist or upper border of the tights (patient specific)
H	Widest part of the pelvis
K*	Length of inner leg, vertically measured from the crotch to the floor (not from measuring point G)
G	2 Finger-widths below the crotch (border must end below the gluteal cleft)
F	Mid-thigh level
E	Middle of the knee-cap (take circumference measurements with knee slightly bent)
D	Head of fibula (2 finger-widths below the knee-cap)
C	At the greatest circumference of the calf
B <sub>1</sub>	Bottom of calf
B	At the slimmest circumference of the lower leg
Y	Measurement across heel and instep with maximum upward flexion of the foot
A	Metatarsophalangeal joint of the little toe
<b>Foot length</b>	Heel to the base of the little toe
<b>Total foot length</b>	Heel to the tip of the longest toe

\* Required information for tights

Buyer Name \_\_\_\_\_


Comments \_\_\_\_\_

Company stamp \_\_\_\_\_

Consignment Patient \_\_\_\_\_

Date \_\_\_\_\_ Customer no. \_\_\_\_\_

Last order no. \_\_\_\_\_

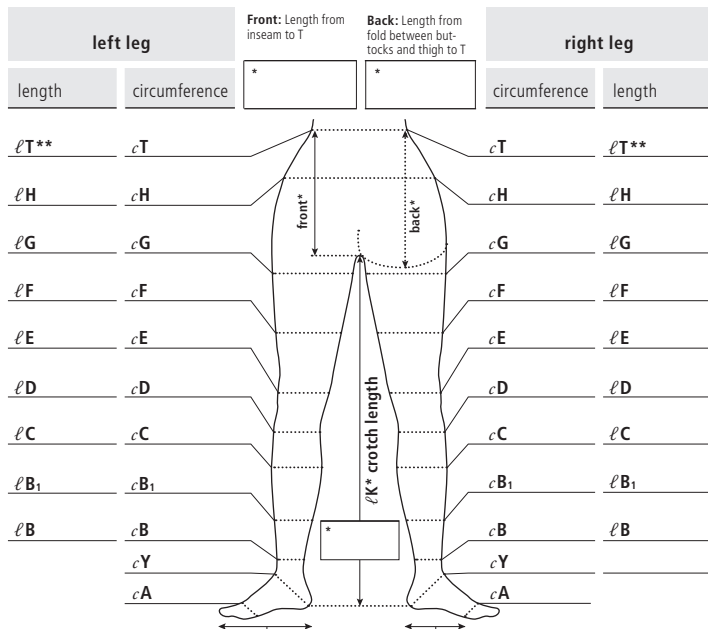

order quantity

piece	pair	whole foot length	shoe size	foot length

- open toe
- closed toe
- asymmetrical closed soft toe<sup>5</sup>

\* Required information for tights  
\*\* Only with AGT (stocking with hip band)

Quality	Ccl.			Colour			
	1	2	3				
memory	<input type="checkbox"/>	<input type="checkbox"/>		<input type="checkbox"/> sugar	<input type="checkbox"/> marzipan <sup>1</sup>	<input type="checkbox"/> honey	<input type="checkbox"/> caramel
				<input type="checkbox"/> brittle	<input type="checkbox"/> nougat <sup>1</sup>	<input type="checkbox"/> chocolate <sup>1</sup>	<input type="checkbox"/> blueberry
				<input type="checkbox"/> poppy seeds <sup>1</sup>	<input type="checkbox"/> liquorice	<input type="checkbox"/> fashion colour: _____ <sup>1</sup>	
memory aloe vera		<input type="checkbox"/>		<input type="checkbox"/> sugar	<input type="checkbox"/> honey	<input type="checkbox"/> caramel	<input type="checkbox"/> brittle
				<input type="checkbox"/> nougat <sup>1</sup>	<input type="checkbox"/> liquorice		
lastofa®	<input type="checkbox"/>	<input type="checkbox"/>		<input type="checkbox"/> mode	<input type="checkbox"/> mode-bright	<input type="checkbox"/> black	<input type="checkbox"/> diamond <sup>4</sup>
				<input type="checkbox"/> nougat <sup>1, 4</sup>	<input type="checkbox"/> ebony <sup>4</sup>		
lastofa® with cotton		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/> mode	<input type="checkbox"/> mode-bright	<input type="checkbox"/> black	

**Grip top**

- lace 5 cm
- silicone 5 cm
- rubber 5 cm
- silicone "Fleur" 5 cm
- silicone 3 cm
- rubber 3 cm
- edge of tights (Trikotrand)<sup>6</sup>

**Design**

- tights AT
- with extra wide panty part ATU
- comfort gusset<sup>2</sup>
- with compressive panty part ATK<sup>2</sup>
- with panty slip ATS<sup>3</sup>
- tights for men ATH<sup>4,6</sup>
- only with AGT (stocking with hip band)<sup>4,6</sup>
- tights single leg ATE<sup>4,6</sup>

<sup>1</sup> 14 days delivery period

<sup>2</sup> Not available with ATS

<sup>3</sup> Not available with Lastofa and Lastofa cotton

<sup>4</sup> Not available with Lastofa Ccl. 1

<sup>5</sup> Not available with Lastofa Cotton

<sup>6</sup> Not available with Memory and Memory Aloe Vera

**Ofa Bamberg GmbH**

Laubanger 20  
D-96052 Bamberg  
Phone + 49 951 6047-777  
export@ofa.de

ofa bamberg